

FINESSE INJECTABLE™

PTA Balloon Dilatation Catheter

SUMMA


Transforming Angioplasty™

Injection ports located 8cm proximal to balloon allow for fluid delivery with the balloon catheter and wire in place.

Angioplasty &
Angiography with
One Device


Ultra-low Crossing Profile

Finesse's extremely low crossing profile and excellent trackability leads to excellent distant lesion crossing and access.

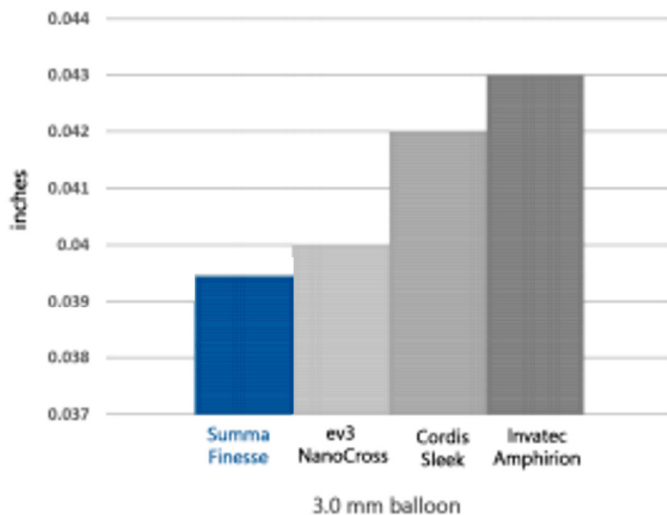

Optimal Lesion Entry Profile

An ultra-low 0.016" tip profile combined with a proprietary hydrophilic coating allows for easy access to the extremities.

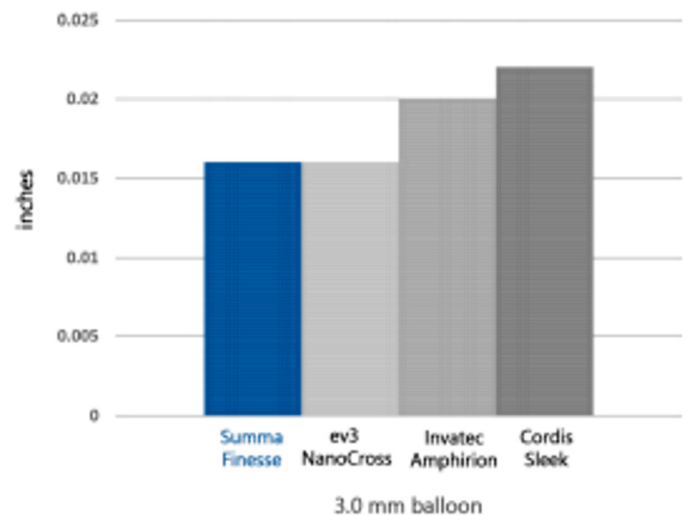

Angioplasty + Angiography. One Device.

The Finesse™ OTW Combination PTA Plus Injection Catheter is the only patented device that combines an 0.014" high pressure balloon catheter with a built-in injection capability that provides targeted injection of contrast agents or therapeutic fluids. This unique design allows physicians to visualize by delivering contrast through the 3 injection ports proximal to the balloon. This can be accomplished while maintaining wire position while the balloon is inflated or deflated.

Technical Specifications

Catheter type	OTW	Balloon lengths	20 – 225 mm
Recommended guide wire	0.014"	Balloon type	Semi-compliant
Tip	Optimal 0.016" entry profile	Distal shaft	2.9F (0.038")
Balloon material	Pebax	Injection ports	3 ports proximal to balloon
Balloon coating	Hydrophilic	Usable length	150 cm
Balloon markers	2 zero profile platinum bands	Introducer sheath compatibility	5F
Balloon diameters	2.0 – 4.0 mm	Guiding catheter compatibility	6F

Ordering Information

Catalog Number	Balloon Diameter (mm)	Balloon Length (mm)	Shaft Length (cm)	Sheath Size (F)	Nominal Pressure (atm)	Rated Burst Pressure (atm)	Guidewire (in)
FS200201508	2	20	150	5	10	12	0.014
FS250201508	2.5	20	150	5	8	12	0.014
FS300201508	3	20	150	5	9	12	0.014
FS350201508	3.5	20	150	5	8	12	0.014
FS400201508	4	20	150	5	7	12	0.014
FS200401508	2	40	150	5	10	12	0.014
FS250401508	2.5	40	150	5	8	12	0.014
FS300401508	3	40	150	5	9	12	0.014
FS350401508	3.5	40	150	5	8	12	0.014
FS400401508	4	40	150	5	7	12	0.014
FS200601508	2	60	150	5	10	12	0.014
FS250601508	2.5	60	150	5	8	12	0.014
FS300601508	3	60	150	5	9	12	0.014
FS350601508	3.5	60	150	5	8	12	0.014
FS400601508	4	60	150	5	7	12	0.014
FS201001508	2	100	150	5	10	12	0.014
FS251001508	2.5	100	150	5	8	12	0.014
FS301001508	3	100	150	5	9	12	0.014
FS351001508	3.5	100	150	5	8	12	0.014
FS401001508	4	100	150	5	7	12	0.014
FS202251508	2	225	150	5	10	12	0.014
FS252251508	2.5	225	150	5	8	12	0.014
FS302251508	3	225	150	5	9	12	0.014
FS352251508	3.5	225	150	5	8	12	0.014
FS402251508	4	225	150	5	7	12	0.014

For Healthcare Professionals only. Important information: Prior to use, refer to the "Instructions For Use" supplied with these devices for indications, contraindications, side effects, suggested procedure, warnings and precautions.

Copyright Summa Therapeutics 2019. All rights reserved. The Summa Therapeutics logo and Finesse logo are trademarks of Summa Therapeutics and may be registered in the U.S. and other countries. All other marks are the property of their respective owners.

Global Headquarters
Summa Therapeutics, LLC
 225 Dyer Street
 Providence, RI 02903
 Tel. 1 (401) 424-1187

www.summaterapeutics.com

SUMMA
 The Specialty Balloon Company